



# RED MAX AJUSTABLE STAND

INSTRUCTION MANUAL

# Table of Contents

<b>Meet Your UltraLux Red Max Adjustable Stand</b>	2
<b>What's in Your Box?</b>	3
<b>How to Assemble Your UltraLux Red Max Adjustable Stand</b>	5
<b>How to Attach the UltraLux Red Max</b>	8
<b>Rotating Your Max Stand</b>	10
<b>Care Instructions</b>	11
<b>Warranty &amp; Contact for Assistance</b>	12

# Meet Your Red Max Adjustable Stand

**Congratulations on amplifying your red light therapy sessions by purchasing the Red Max Adjustable Stand!**

**The adjustable stand is designed to enhance your experience with the UltraLux Red Max device, allowing you to use it in your preferred comfort zone.**

Not only is the Max Stand made of durable metal, but it includes wheels for increased mobility. Unlike other red light therapy stands, the Red Max Adjustable Stand offers both vertical and horizontal tilt options, providing true tailored red light sessions to suit your preferences.



# What's in Your Box?



8x **A**



1x **B**



4x **C**



1x **D**



2x **E**



16x **F**



6x **G**




4x **H**





1x **I**



4x **J**
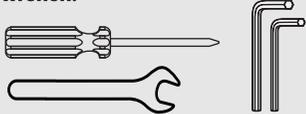
<b>A</b>	Base and leg bolts	<b>H</b>	Washers	<b>O</b>	Panel mounting brackets
<b>B</b>	Height Adjusting Hand Knob	<b>I</b>	Lynch pin	<p>Every part of the Max Stand is labeled with a corresponding letter of the alphabet.</p> <p><b>All parts can be assembled using the included screwdriver and wrench.</b></p> 	
<b>C</b>	Castors	<b>J</b>	Base bolts covers		
<b>D</b>	L-Bar hand knob	<b>K</b>	Base T-Bar		
<b>E</b>	L-Bar bolts	<b>L</b>	Bracket swing arm		
<b>F</b>	Panel mounting screws	<b>M</b>	Base legs		
<b>G</b>	Bracket mounting screws	<b>N</b>	L-Bar		

# What's in Your Box?

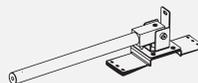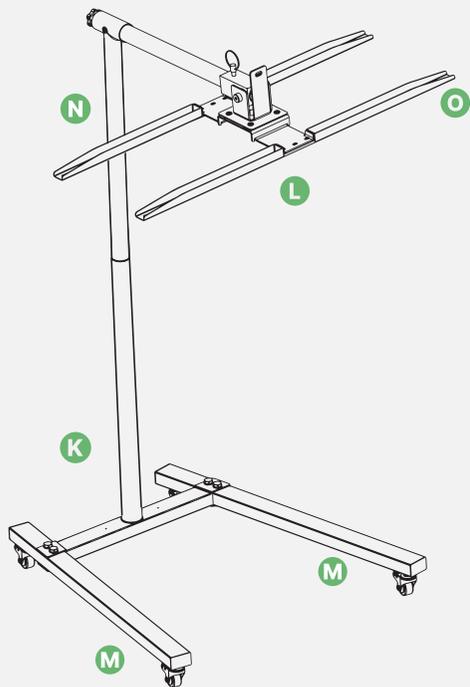


# How to Assemble



FOR VIDEO INSTRUCTIONS  
SCAN THIS QR CODE WITH  
YOUR PHONE'S CAMERA


1. Begin by making sure you have the Base T-Bar (K), Base Legs (M), Base and Leg Bolts (A), Base Bolt Covers (J), and Washers (H). Match the holes of the T-Bar (K), to the holes of the Base Legs (M).

Use the 4 Base and Leg Bolts (A) and 4 Washers (H) to screw the T-Bar (K) onto the Base Legs (M). Then place the Base Bolt Covers (J) onto the screwed Bolts.

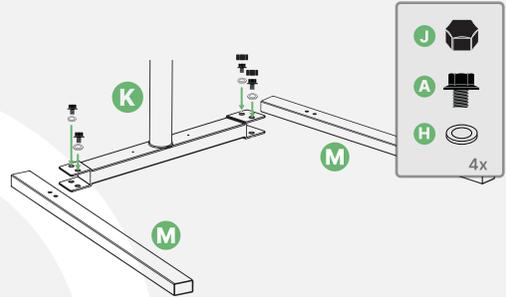

2. Next screw the remaining Base and Leg Bolts (A) to the Base Legs (M).


Then screw in the Castors (C) onto the ends of both Base-Legs (M).

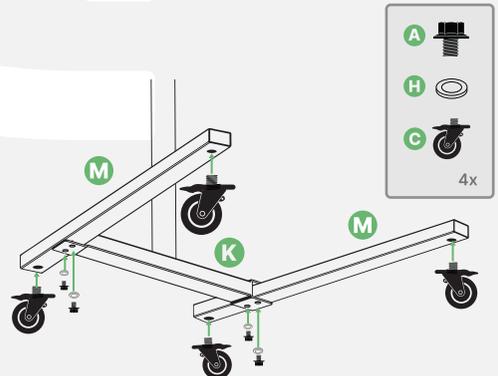


# How to Assemble

**3.** Grab the L-Bar (N) and slide the end of Bracket Swing Arm (L) into the opening end of the L-Bar (N). Screw the Bracket Swing Arm (L) in place, using the L-Bar Bolts (E).

Next, screw in the L-Bar hand knob (D) to the back of the L-Bar (N) to ensure a tight fit. Then place the Lynch Pin (I) to the back of the Bracket Swing Arm (L) to lock in place.

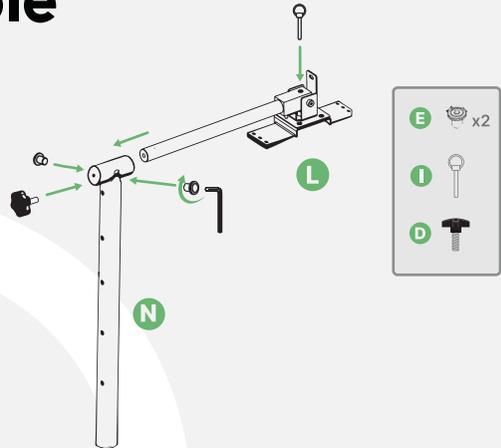


## Adjusting the Position Red Max

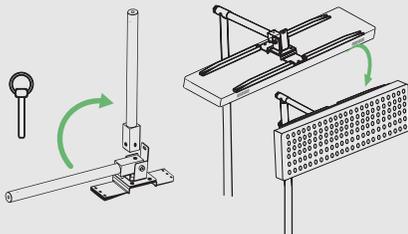


Place your hand under your Red Max to slightly lift "up" and lessen the tension on the Lynch Pin (I) to be able to remove it. Doing so will allow you to tilt the Bracket Swing Arm (L) up or down. When you have your desired position, reinsert the Lynch Pin (I) to lock it in place.

# How to Assemble

- 4.** Then use the 6 Bracket Mounting Screws (G) to attach the Panel Mounting Brackets (O) to both ends of the Bracket Swing Arm (L), lining up each hole with the holes in the middle of the Panel Mounting Brackets (O).

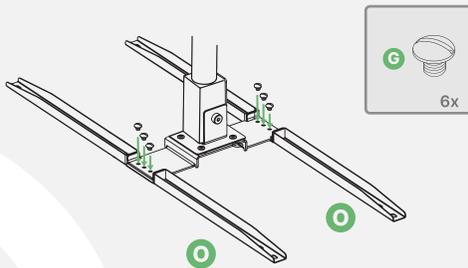


- 5.** Finally, attach the L-Bar (N) to the opening of the T-Bar (K) screwing in the Height Adjusting hand knob (B) into the lined up holes.

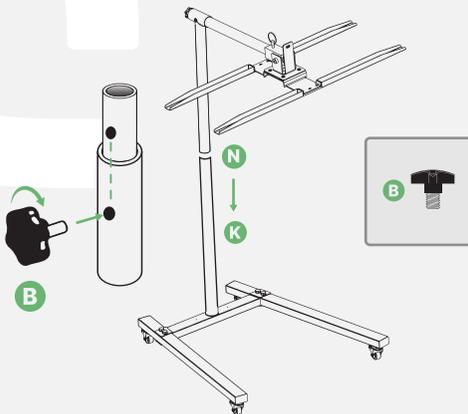


# How to Attach the Red Max

1. Before attaching the Red Max to the Stand, we recommend removing the rubber feet from the back of the Red Max.

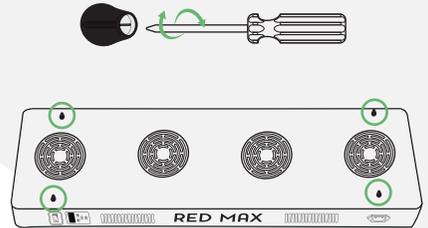



2. To avoid damage, remove the L-Bar (N) from the Base T-Bar (K) before attaching the Red Max.

Then use the 16 Panel Mounting Screws (F) to screw the back of the UltraLux Red Max panel to the Panel Mounting Brackets (O).

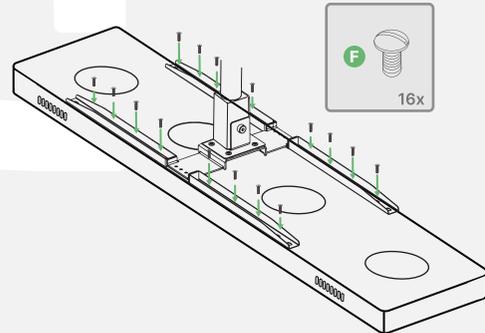

# How to Attach the Red Max


- 3.** Lastly, reattach the L-Bar (N) to the Base T-Bar (K), using the Height Adjusting Hand knob (B), and your Red Max is ready to use!

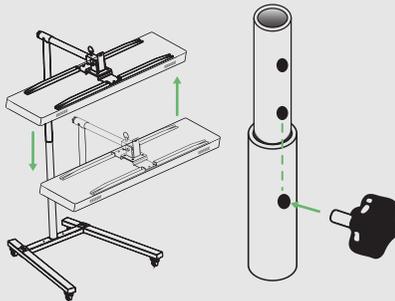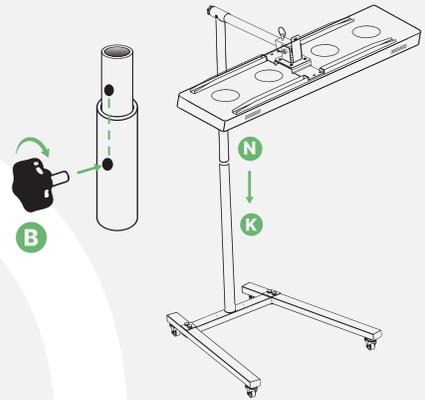

## Stand Height Adjustment

While supporting the weight of the L-Bar (N), unscrew the Height Adjusting Hand Knob (B) from the Base T-Bar (K). You can then lower or raise the L-Bar (N) to the desired height where 2 holes align on the Base T-Bar (K). When you have the desired position, screw the Height Adjusting Hand Knob (B) back into the Base T-Bar (K).


## Rotating Your Max Stand

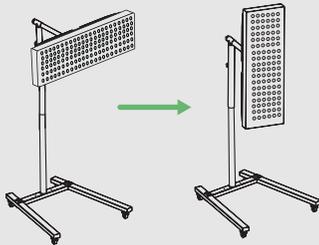

You can rotate your Max Stand to be positioned vertically or horizontally. All instructions above were shown in the horizontal position with the ability to tilt the unit to be used while laying down or facing outwards.

You may rotate your Max Stand 90° without needing to remove your Red Max Panel from the Brackets. Simply unscrew the 4 screws on the Bracket Swing Arm (L) and rotate the entire section holding the Red Max.

**Do not worry - the device will not fall when you unscrew the 4 screws.**

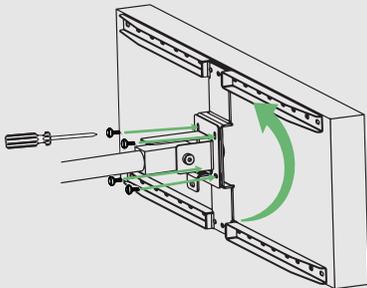

There is a large pre-installed bolt holding the pieces together on a swivel. Once you have rotated the Bracket 90°, simply reposition the white plastic spacers with a screw 1 by 1 and tighten them.

# Care Instructions

When storing the Max Stand with the Red Max attached, it is important to make sure the Red Max is unplugged, before keeping it in a cool, dry area away from direct sunlight.

If the Red Max is not attached to the Max Stand, feel free to store your Max Stand wherever is convenient, without compromising its structure.

## To Clean the Red Max Stand:

Use a damp cloth or paper towel to wipe the Max Stand.

Ensure the stand is completely dry before attaching the Red Max.



# Warranty

The UltraLux Red Max Adjustable Stand comes with a **1-year manufacturer warranty, and only valid for the original buyer.** In order for the warranty to apply, the item must be in like new condition, excluding the defect.

Contact **support@ultraluxhealth.com** with photo or video evidence of the defective item. Whether or not the item is covered under warranty is under full discretion of UltraLux Health. If the item is covered, a replacement item will be shipped to the customer at no cost to the customer. Warranties cover labor, defective parts, and the return shipping of the new/repaired product to the customer.

# Contact For Assistance

Now that you know exactly how to operate and use your UltraLux Red Max Adjustable Stand, it's time to take full advantage of the amazing benefits!

**Feel free to reach out to our customer service team with any questions or concerns**



✉ Support@ultraluxhealth.com    💬 Live Chat : ultraluxhealth.com

☎ 385-360-9755 or Internationally Toll Free at: 855-479-4149

M-F 9:00am - 5:00pm MST



**Questions? We are happy to help.**

✉ [Support@ultraluxhealth.com](mailto:Support@ultraluxhealth.com)  Live Chat : [ultraluxhealth.com](https://ultraluxhealth.com)

 385-360-9755 or Internationally Toll Free at: 855-479-4149

M-F 9:00am - 5:00pm MST