



RED PRO

INSTRUCTION MANUAL




Table of Contents

Meet Your UltraLux Red Pro	2
Health Benefits from Using the Red Pro	4
What's In Your Red Pro Box?	5
How to Use Your Red Pro	6
Safety Precautions & Care Instructions	7
Frequently Asked Questions (FAQ)	9
Warranty & Contact for Assistance	11
Red Pro User Health Guide	13
How To Optimize Red Pro Benefits	15
What To Avoid	16

Meet Your UltraLux Red Pro

Congratulations on taking this next step in your health journey! We hope you'll enjoy the expert-level restorative and healing aspects of red light therapy in your in your own home or on the go.

While red light therapy is renowned for its skin benefits, **the UltraLux Red Pro harnesses red light wavelengths proven to increase muscle growth and provide relief from pain and inflammation.** With the power of natural healing in your palm, achieving your healthiest self has never been easier.

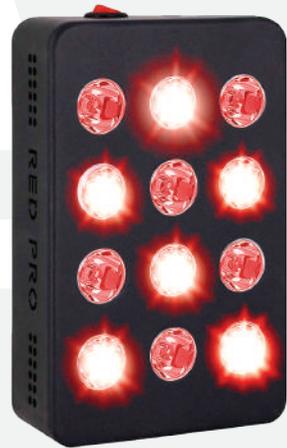

What Sets the UltraLux Red Pro Apart?

Powerful Single Chip LEDs: Our custom LEDs provide the highest irradiance output for any handheld device because of their single chip design. This allows for one bulb to maximize its invisible infrared or red light.

Clinic-Quality Healing at Home

Experience the best of healing technology without the hassle of setting up a doctor's appointment. Rejuvenate your skin and support your body's healing without leaving the comfort of your home.

FDA-Cleared Wavelengths for Weight Loss

The UltraLux Red Pro includes red and infrared light wavelengths to heal the effects of aging and skin conditions.



UltraLux Red Pro Benefits You May Experience

The UltraLux Red Wrap combines red (660nm) and infrared (850nm) wavelengths which are ideal for providing therapeutic relief.

Nanometers*



630nm


Anti-aging
Wound and injury healing
Pain relief
Scar and skin condition healing
Stem cell regeneration

850nm

Tightens skin
Muscle building
Increased athletic performance
Post-workout recovery

*Nanometers are the unit of measure for light wavelengths.

What's in Your Red Pro Box?

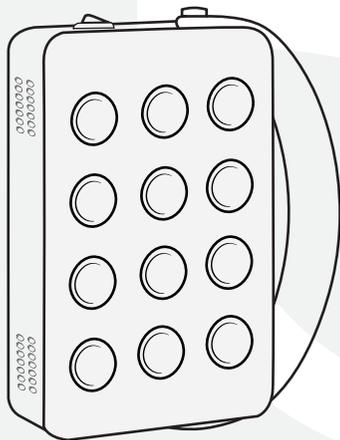

1x UltraLux Red Pro
With 2 Silicon
Hand Straps



1x USB
Cable




1x Adaptor



1x Red Pro
Stand



1x Light Protection
Goggles

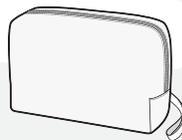


1x Hard Cover
Carry Case




1x User
Manual


How to Use Your Red Pro

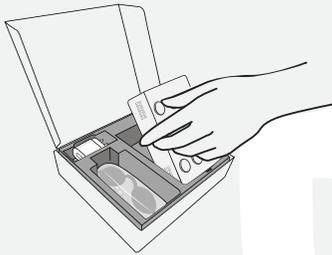


1. Take out the Red Pro and all contents from the box.

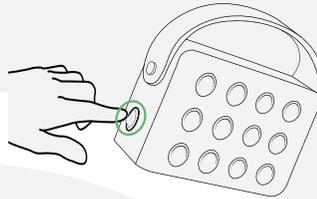

2. Ensure your Red Pro is charged before use and turn on the power switch located on the side of the Red Pro.

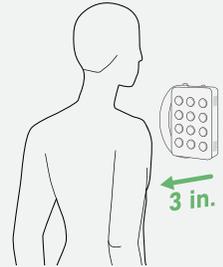

3. Hover the Red Pro approximately 3 inches away from your desired treatment area.

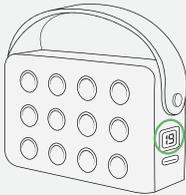


4. The digital screen on the side of the device will automatically set a timer for 20 minutes




**FOR VIDEO INSTRUCTIONS
SCAN THIS QR CODE WITH
YOUR PHONE'S CAMERA**

Charging Your Red Pro

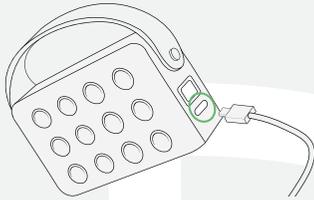


1. Plug the provided USB-C charging cord into the charging port on the side of the device, and the USB end of the charging cord into the wall adaptor.

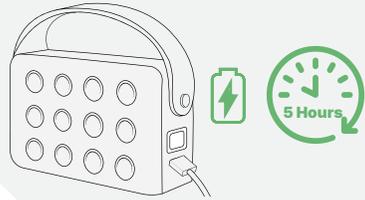


2. Charge your device for a few hours to ensure a full charge.

Note: Do not turn on the lights when the device is charging. It will decrease the battery lifespan.

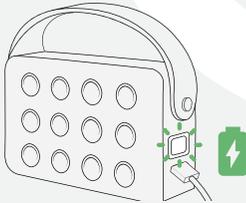

3. You will know the device is fully charged when the digital timer interface goes from flashing to continuously lit.

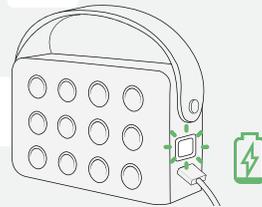

4. A flashing interface will also indicate the device is low battery and needs to be charged.

Safety Instructions



This device is water-resistant, **NOT WATERPROOF.**

Do not submerge the Red Pro or its components in water.

- Not intended for use by children or pregnant women.
- Do not exceed 20-minute sessions.
- Do not use if damaged as it could result in electric shock.
- Do not use in conjunction with sunscreen, ointments, or balms.
- Do not use on or near hot surfaces.
- Do not modify or attempt to modify this device.

Consult Your Doctor if You:

- Are pregnant or breastfeeding.
- Recently had a steroid injection or use topical steroids.
- Have any medical questions regarding red light therapy.
- Experience negative reactions or effects from using this device.
- Are taking drugs or medication known to cause light sensitivity.

Care Instructions

Before storing your Red Pro, unplug it and keep it in the carrying case in a cool, dry area away from direct sunlight.

To clean the Red Pro:

- Make sure all cords are disconnected, and that the panel is powered off.
- Use a damp cloth or paper towel to clean the panel. Only clean the outside portion of the Red Pro, avoiding the LED diodes.
- Never submerge your panel in water or other liquids.
- Allow your Red Pro to air dry before using it again.

The statements have not been evaluated by the FDA (U.S. Food & Drug Administration). Our products are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness / fitness products that do not require FDA clearance, in accordance with the "General Wellness: Policy on Low Risk Devices" draft released January 20, 2015

Frequently Asked Questions

How does red light therapy work?

Red light therapy delivers natural, non-invasive wavelengths of light to the cells in our bodies. When these light wavelengths are absorbed, the mitochondrion are stimulated, creating more energy within our bodies. As a result, this cellular energy is used to increase the health and effectiveness of our bodily functions.

Are there side effects from using this device?

Red light therapy is not associated with any adverse side effects if used as directed.

How long will it take to notice results?

You may start noticing visible or physical results in as little as 4 weeks of consistent use. Individual results may vary.

Is Red Light bad for my eyes?

Red light therapy is shown to improve eye health! Most customers find it useful to close their eyes while using the Red Pro.

What is irradiance?

Irradiance measures the strength and effectiveness of a light wavelength. UltraLux Health uses high irradiances on all red light therapy products to ensure real, powerful results.

Will I get a tan from the red light?

Irradiance measures the strength and effectiveness of a light wavelength. UltraLux Health uses high irradiances on all red light therapy products to ensure real, powerful results.

Frequently Asked Questions

Do I need to clean the device? How do I clean it?

Red light therapy delivers natural, non-invasive wavelengths of light to the cells in our bodies. When these light wavelengths are absorbed, the mitochondrion are stimulated, creating more energy within our bodies. As a result, this cellular energy is used to increase the health and effectiveness of our bodily functions.

Does red light therapy hurt?

No, it should not hurt. Red light therapy should feel like a comfortable, slight warmth. If you experience skin irritation or pain, stop using and consult a medical professional right away.

Do I need a prescription from my doctor to use the Red Pro?

Not at all! Our light therapy device uses FDA-cleared wavelengths that can be bought over the counter for anyone experiencing pain issues. If you have any health concerns, or light sensitivity issues, please consult a medical professional before use.

Is the Red Pro battery operated?

The Red Pro includes a powerful, long-lasting rechargeable battery, so no need to purchase batteries separately! The Red Pro requires a full charge before use, but after that, it is completely mobile for your treatment session.

How often can I use the Red Pro?

For best results, use the Red Pro for 1-2, 15-minute sessions per day for a minimum of 4 weeks. We strongly advise you to only use the Red Pro for these recommended time allotments. Do not exceed 30 minutes per treatment session.

Warranty

The UltraLux Red Pro comes with a **2-year manufacturer warranty, and only valid for the original buyer.** In order for the warranty to apply, the item must be in like new condition, excluding the defect.

Contact **support@ultraluxhealth.com** with photo or video evidence of the defective item. Whether or not the item is covered under warranty is under full discretion of UltraLux Health. If the item is covered, a replacement item will be shipped to the customer at no cost to the customer. Warranties cover labor, defective parts, and the return shipping of the new/repaired product to the customer.

Contact For Assistance

Now that you know exactly how to operate and use your UltraLux Red Cap, it's time to take full advantage of the amazing benefits!

Feel free to reach out to our customer service team with any questions or concerns



✉ Support@ultraluxhealth.com 💬 Live Chat : ultraluxhealth.com

☎ 385-360-9755 or Internationally Toll Free at: 855-479-4149

M-F 9:00am - 5:00pm MST



RED PRO

HEALTH GUIDE



UltraLux Red Pro Health Guide

The UltraLux Health Red Pro has benefits for all types of skin, muscle, and hair conditions. Whether desiring to achieve healthier skin or simply want to improve muscle recovery, the Red Pro can be your solution. This Red Pro Health Guide will provide you with information to facilitate quality results.

Supported Principles of the Red Pro

Red light therapy (RLT) supports various principles related to internal biological processes that promote body, skin, and mental health. The UltraLux Red Pro supports these same processes while targeting the skin and muscles. Additionally, we advise consulting a doctor if you have questions about a specific condition related to these targeted areas. However, if your specific condition would improve by experiencing a principle listed below, then the Red Pro is right for you!



Increased Collagen

Collagen is a structural protein that provides support to the skin and hair. Red light therapy's ability to boost collagen production may contribute to accelerated healing of wounds, scars, wrinkles, and skin injuries.



Stimulation of Cellular Activity

Red light therapy enhances cellular function by promoting increased energy production (ATP) within the cells. This may lead to improved cell metabolism and overall cellular health in the body.



Reduction of Inflammation

Red light therapy enhances cellular function by promoting increased energy production (ATP) within the cells. This may lead to improved cell metabolism and overall cellular health in the body.



Increased Blood Circulation

Red light therapy stimulates blood flow which aids in the removal of waste products and toxins from tissues. Improved blood circulation ensures a more efficient delivery of oxygen, nutrients, and other essential substances to the cells.



Mental Health

Exposure to red light is believed to stimulate cellular activity, particularly in the mitochondria, leading to enhanced energy production and improved neurotransmitter function. The release of certain neurotransmitters, such as serotonin and endorphins, which play key roles in mood regulation and emotional well-being.

How to Optimize Your Red Pro Benefits

Tailoring your red light therapy sessions is great for achieving goals and maximizing the potential benefits of this treatment faster. For the best results, we recommend the following:



Ensure skin is free from makeup, lotions, sunscreens, or other barriers.

Clean the target body area with a gentle cleanser before your red light therapy session.



Remove Clothing or other barriers from the target body area.

Bare skin absorbs the red light better than covered skin.



What to Avoid

When using the UltraLux Red Pro, it is important to maintain your health and safety. We recommend avoiding the following in order to receive the results you want, while protecting yourself.

Overexposure

We recommend not exceeding 30 minutes per session. Do not use your Red Pro overnight.

Not Consulting a Doctor About Conditions and Medications

If you suffer from a specific scalp condition or are taking photosensitizing medications, and do not consult with a healthcare professional or dermatologist, there may be adverse effects. Talk to your doctor, and they can provide personalized advice based on your health status.

Inconsistency in Use

Consistency is key with red light therapy. Avoid sporadic or irregular usage, as regular and consistent sessions are more effective for achieving quality results. But do not exceed recommended session time.

DISCLAIMER : The information provided in this health guide is intended for informational purposes only and is not a substitute for professional medical advice, diagnosis, or medical treatment. Use the information in this guide responsibly, and seek professional medical guidance from your doctor if you have questions concerning your specific health concerns.

Product Specifications

Product Name	UltraLux Red Pro Red Light Therapy Panel
Irradiance	6 inches 93n mw/cm2 4 inches 173 mw/cm2
Wavelength	660nm / 850nm
Charge Time	Approximately 5 Hours
Sessions per Charge	2-3
Size	6 x 3.7 x 1.5 inches
Timer	20 minutes
Voltage	DC 5V 2A
LED Count	6 Single Chip NIR LEDs 6 Single Chip RED



Questions? We are happy to help.

✉ Support@ultraluxhealth.com  Live Chat : ultraluxhealth.com

☎ 385-360-9755 or Internationally Toll Free at: 855-479-4149

M-F 9:00am - 5:00pm MST